


STACKED BIOCASSETTE

for cassette printers



CODE	COLOUR	PACKAGING
07-9700	White	1.600 pcs (40 stacks of 40 pcs – taped)
07-9710	Orange	1.600 pcs (40 stacks of 40 pcs – taped)
07-9720	Blue	1.600 pcs (40 stacks of 40 pcs – taped)
07-9730	Yellow	1.600 pcs (40 stacks of 40 pcs – taped)
07-9740	Lilac	1.600 pcs (40 stacks of 40 pcs – taped)
07-9750	Pink	1.600 pcs (40 stacks of 40 pcs – taped)
07-9760	Green	1.600 pcs (40 stacks of 40 pcs – taped)

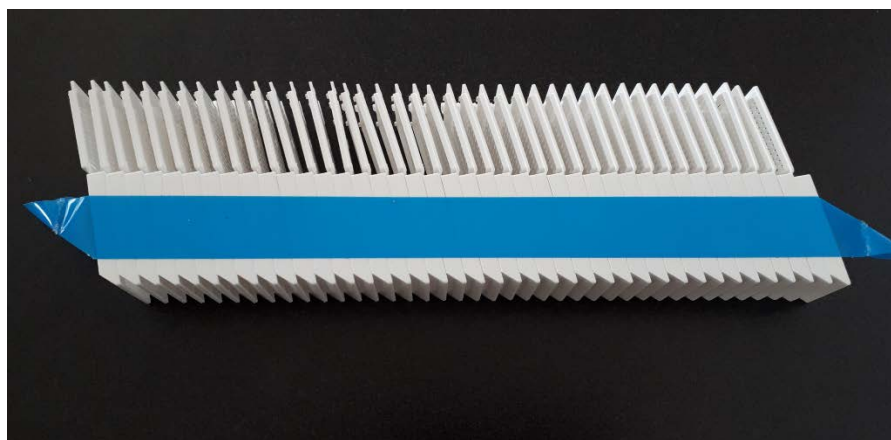

In vitro diagnostic – medical device



Manufacturer: Bio-Optica Milano S.p.A.



Disposable product

Date of issue: 10/01/2019

Rev. 03

TECHNICAL FEATURES

Standard Embedding cassettes with double closure system to guarantee a perfect resistance during tissue processing. Suitable surface for automatic cassette printers. The cassettes are stacked with tape, with attached lid.





Material	Acetalic polymer (polyoximethylene polyacetal). The acetalic polymer is stable to the high temperatures, to microwave operations, to solvents, to acids and alkalis; then this material is optimal for the processing and the embedding passages in the histological process.
Internal size	30 x 25 x 5 mm
External size	40 x 28 x 6 mm
Holes size	(5.8 ± 0.2) x (1 ± 0.1) mm
Surface inclination	45°
Medium weight for each cassette	3,40 g
Internal packaging	Carton box
External packaging	Carton box
Expiry date	Not applicable for this product
Storage precautions	Not necessary for this kind of product
Warning	Not provided for this kind of product
Disposal	Dispose in conformity with local law

The sample's identification is possible thanks to the smooth surface on the frontal cassette's side. This surface is fitted to the writing with automatic cassette printers. Not suitable for handwriting (pencil or pen).

The embedding cassettes are compatible with the most common filing system on the market.

Bio cassettes are manufactured with our best technical skills, therefore, due to the complexity of the details and the type of material used, is to be considered a physiological quantity of defective cassettes equal to 5,000 PPM (Parts Per Million, percentage of 0.5%). Within this amount, Bio-Optica does not accept complaints and / or requests for warranty replacement.

COMPATIBILITY WITH THE WRITING SYSTEMS ON THE MARKET

Code	 Primera	 Leica IPC	 Sakura SmartWrite	 Sakura AutoWrite	 Thermo PrintMate	 Hand writing
07-9700	✓	✓	✓	✓	✗	✗
07-9710	✓	✓	✓	✓	✗	✗
07-9720	✓	✓	✓	✓	✗	✗
07-9730	✓	✓	✓	✓	✗	✗
07-9740	✓	✓	✓	✓	✗	✗
07-9750	✓	✓	✓	✓	✗	✗
07-9760	✓	✓	✓	✓	✗	✗

Date of issue: 10/01/2019

Rev. 03