

OSOM[®]

Mono Test

TO AID IN THE DETECTION OF MONONUCLEOSIS



FOR USE WITH WHOLE BLOOD, SERUM, OR PLASMA

FEATURES:



- ✓ **ACCURATE;** >99% sensitivity; 96% specificity
- ✓ **FAST;** Results in 5 minutes
- ✓ **NO AGE LIMITATIONS**
- ✓ **EASY-TO-READ;** Two-color results
- ✓ **WORLD-CLASS TECHNICAL SUPPORT;** Staffed with certified medical technologists and skilled professionals



TEST PROCEDURE / RESULTS



For serum, plasma or whole blood from collection tubes, add 1 drop to test tube from pipette.

For fingertip blood, fill capillary tube and dispense blood into test tube.




Add 1 drop of diluent to test tube. Mix.

Place the test stick into the treated sample.



Read results at 5 MINUTES.

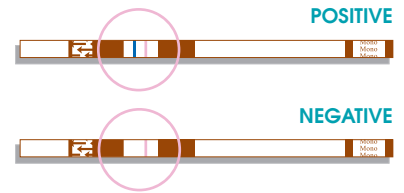

ORDERING INFORMATION

○ OSOM[®] MONO TEST

CONFIGURATION

25 TESTS

CATALOG NUMBER

145 OR 145E (CE MARKED)



Mediq Suomi Oy
Luomanportti 3, PL 115
02201 Espoo, puh. 020 112 1500
www.mediq.fi



INTERNATIONAL
Sekisui Diagnostics (UK) Limited
50 Gibson Drive
Kings Hill, West Malling, Kent, ME19 4AF, UK
Phone: +44 (0) 1732 220022
Fax: +44 (0) 1732 220024
Email: info@sekisuidiagnostics.com
ISO 9001:2000 Certified, ISO 13485:2003 Certified

