




Versatus™-W

Intravenous catheter with wings /
and injection port

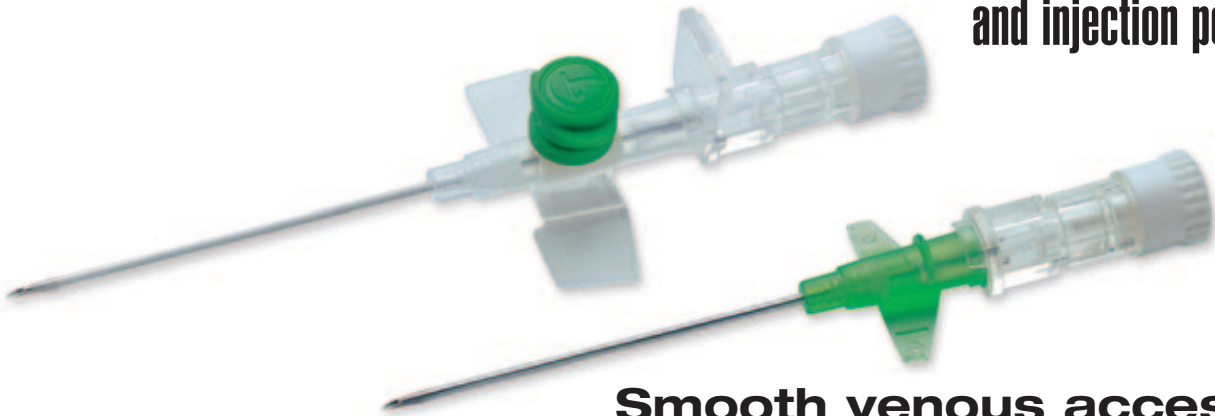

**Smooth venous access
Better comfort for the patient**

Material

The catheter is made from polyurethane (PUR) which is firm at room temperature and softer at body temperature, kink resistant and biocompatible: patient comfort is enhanced. This catheter is radio-opaque.

Internal needle

The quality of the bevel of the internal needle ensures a precise incision.

Conical catheter tip

The perfectly tapered catheter tip allows the catheter to enter the vein without causing any venous trauma and allows for an easy placement. In this way the comfort of the patient is guaranteed. The optimal bevel/catheter tip distance reduces the risk of passing through the vein.





Winged catheter hub

The wings ensure easy fixation of the catheter. The EN/ISO coloured hub of the VERSATUS™-W with wings confirms the catheter size.

Injection port

This allows the intermittent injection of medicines during infusion. Due to central position, there are less risks for damaging the vein and patient trauma during intermittent injection.

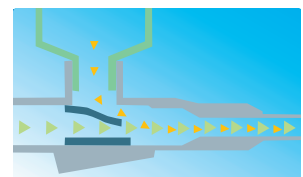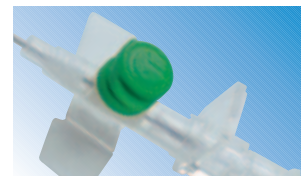

Flashback

The flashback chamber is transparent for easy detection of the flashback. A venting plug, having an hydrophobic filter, allows the passage of air, ensures reliable closure to avoid any risk of



contamination. The Luer Lock cap allows the extremity of the catheter to be closed after withdrawal of the internal needle.

Obturator

A Luer Lock obturator with the same colour code as the catheter is available for each gauge (except gauge 26). This prevents the formation of clots in the lumen of the catheter due to any interrupted use.



VERSATUS™-W with wings and injection port References	EN/ISO colour code	External catheter diameter		Internal catheter diameter mm	Catheter length mm	Flow rate ml/min	Obturator
		Gauge	mm				
SR+DU2619PX	Violet	26	0.64	0.45	19	19	-----
SR+DU2419PX	Yellow	24	0.74	0.55	19	29	SO+DM2419P
SR+DU2225PX	Blue	22	0.90	0.65	25	42	SO+DM2225P
SR+DU2032PX	Pink	20	1.00	0.75	32	59	SO+DM2032P
SR+DU1832PX	Green	18	1.30	0.95	32	103	SO+DM1832P
SR+DU1845PX	Green	18	1.30	0.95	45	96	SO+DM1845P
SR+DU1745PX	White	17	1.50	1.15	45	155	SO+DM1745P
SR+DU1645PX	Grey	16	1.75	1.35	45	225	SO+DM1645P
SR+DU1445PX	Orange	14	2.00	1.55	45	290	SO+DM1445P

VERSATUS™-W with wings References	EN/ISO colour code	External catheter diameter		Internal catheter diameter mm	Catheter length mm	Flow rate ml/min	Obturator
		Gauge	mm				
SR+DU2619WX	Violet	26	0.64	0.45	19	19	-----
SR+DU2419WX	Yellow	24	0.74	0.55	19	29	SO+DM2419W
SR+DU2225WX	Blue	22	0.90	0.65	25	42	SO+DM2225W
SR+DU2032WX	Pink	20	1.00	0.75	32	59	SO+DM2032W
SR+DU1832WX	Green	18	1.30	0.95	32	103	SO+DM1832W
SR+DU1845WX	Green	18	1.30	0.95	45	96	SO+DM1845W
SR+DU1745WX	White	17	1.50	1.15	45	155	SO+DM1745W
SR+DU1645WX	Grey	16	1.75	1.35	45	225	SO+DM1645W
SR+DU1445WX	Orange	14	2.00	1.55	45	290	SO+DM1445W



Every day,
as healthcare
professionals,
you are working
to ensure life
and to contribute to
the well-being
of your patients.

Every day,
Terumo Europe
provides you with
the clinical systems,
means and solutions
to keep life flowing.



We keep life flowing

TERUMO®

TERUMO CORPORATION
44-1, 2-chome
Hatahaya
Shibuya-ku
Tokyo 151-0072
JAPAN
Tel. +81/3 3374.8111
Fax +81/3 3374.8196

TERUMO EUROPE N.V.
Researchpark Zone 2 Haasrode
Interleuvenlaan 40
B-3001 Leuven
BELGIUM
Tel. +32/16.38.12.11
Fax +32/16.40.02.49

TERUMO EUROPE N.V.
BENELUX SALES DIVISION
Researchpark Zone 2 Haasrode
Interleuvenlaan 40
B-3001 Leuven
BELGIUM
Tel. +32/16.39.25.80
Fax +32/16.39.25.99
THE NETHERLANDS
Tel. 0800/022.03.96
Fax 0800/022.04.14

LABORATOIRES TERUMO
FRANCE S.A.
1, Parc Ariane
Bâtiment Uranus
Rue Hélène Boucher
78284 Guyancourt Cedex
FRANCE
Tel. +33/1 30.96.13.00
Fax +33/1 30.43.60.85

TERUMO DEUTSCHLAND GmbH
Hauptstrasse 87
65760 Eschborn
GERMANY
Tel. +49/61.96.80.230
Fax +49/61.96.80.23.200

TERUMO EUROPE ESPAÑA SL
Avenida de Burgos 16-D
Edificio EUROMOR
Madrid 28036
ESPAÑA
Tel. +34/91 383.90.69
+34/91 383.91.61
Fax +34/91 767.04.60

TERUMO ITALIA SRL
Via Simone Martini 143/145
I-00142 Roma
ITALIA
Tel. +39/06.519.61.420
Fax +39/06.503.04.07

TERUMO SWEDEN AB
Nya Varvet
Byggnad 90
SE-426 71 Västra Frölunda
SWEDEN
Tel. +46/31.748.58.80
Fax +46/31.748.58.90

TERUMO DENMARK
Filial af TERUMO SWEDEN AB
Tel. +45/7020.93.80
Fax +45/7020.94.80

TERUMO UK Ltd.
Tamesis
The Causeway
Egham
Surrey
TW20 9AW
UNITED KINGDOM
Tel. +44/1784.476.200
Fax +44/1784.476.201

®: Registered Trademark
Published by Terumo Europe N.V.
www.terumo.com
www.terumo-europe.com